





DIOLASHeer ice -1200Pro

Latest Double Handles Diode Laser Hair Removal



- Comfortable
- Safety
- Effective

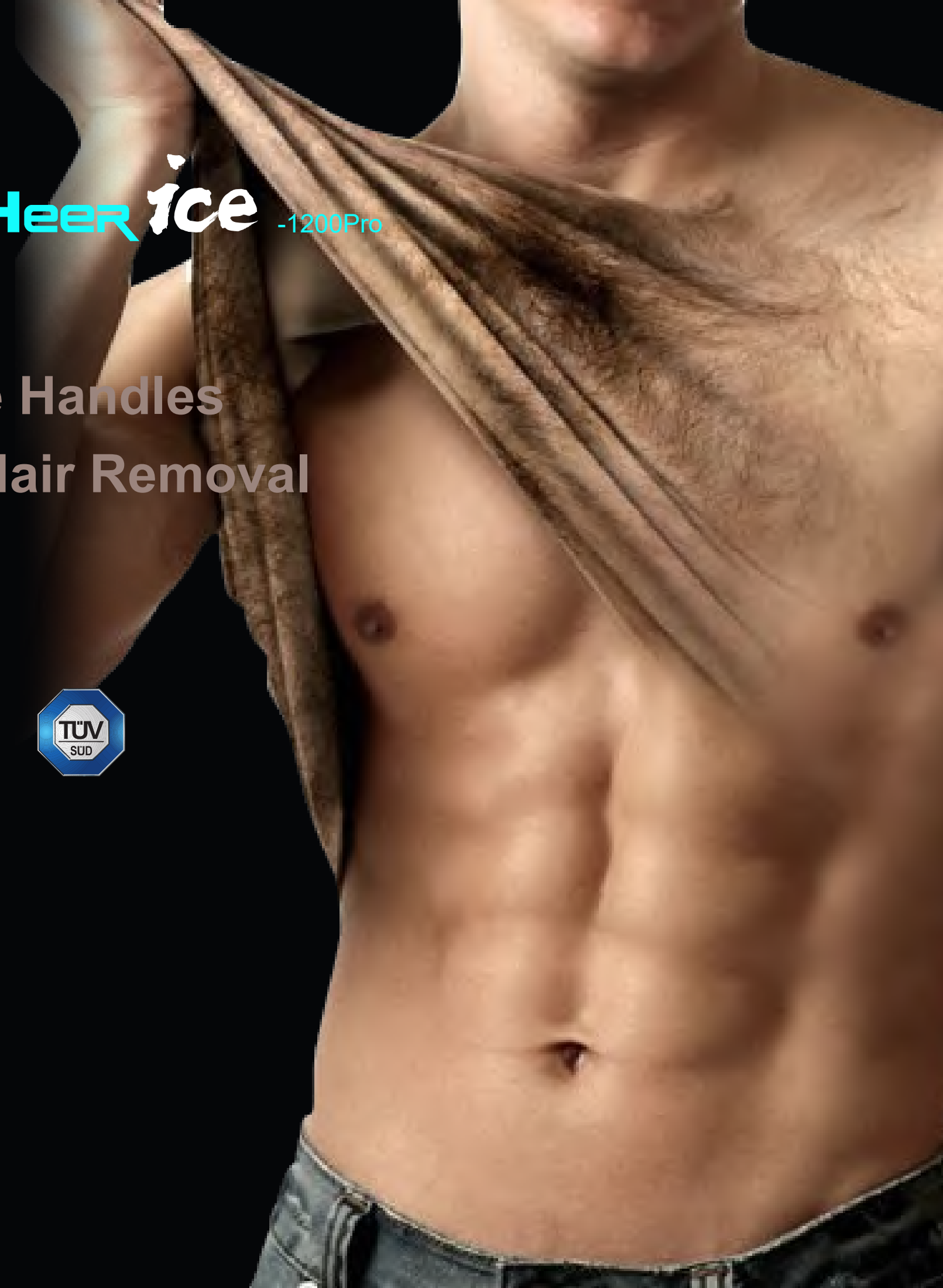

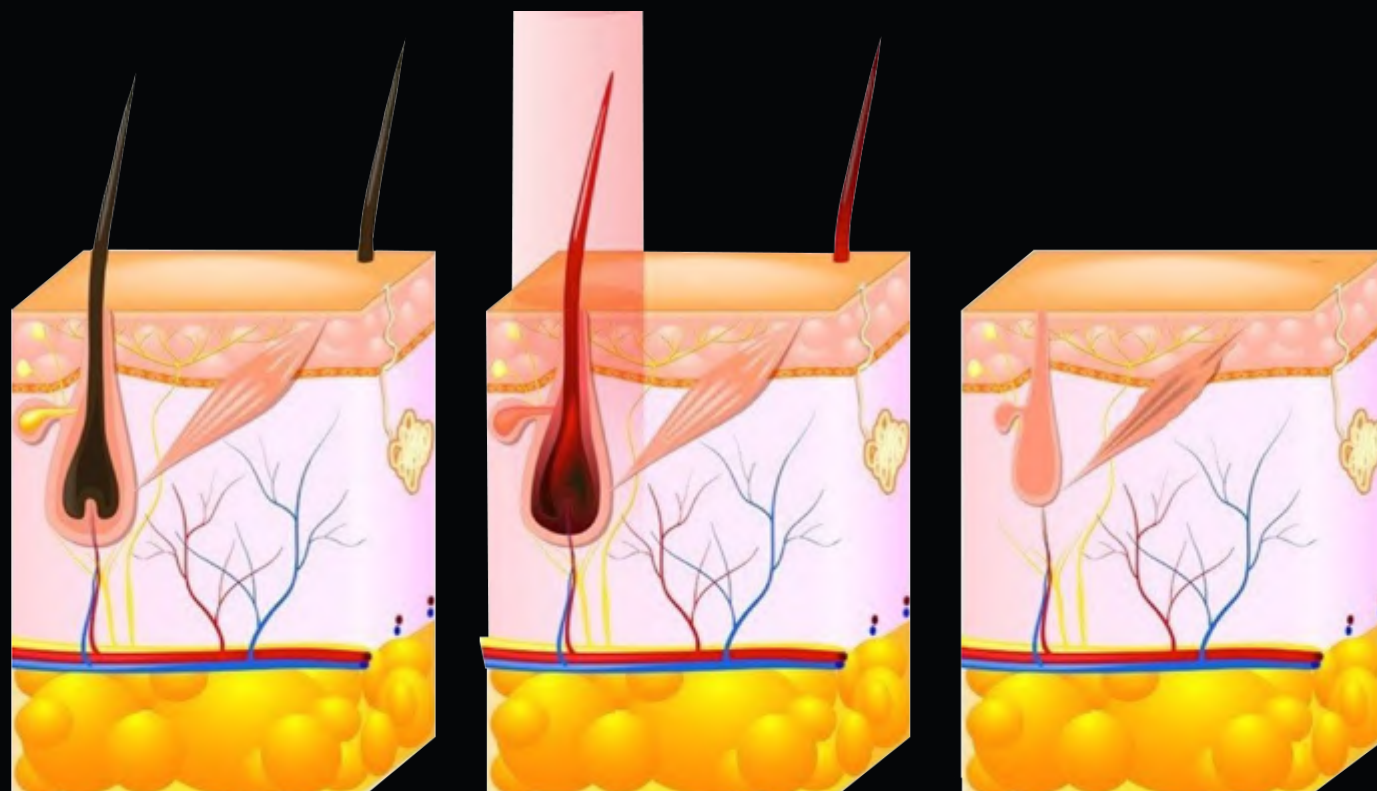


Before

The melanin absorbs
the laser and heats up

After

PRINCIPLE OF TREATMENT

The principle of diode laser hair removal treatment is based on selective photothermal decomposition. The melanosomes in the hair follicle can selectively absorb the energy of the laser. The laser energy emitted by the machine is easily absorbed by the colored hair follicle without damaging the epidermal tissue. The energy emitted by the laser is absorbed by the hair and the pigment in the hair follicle. It is converted into heat, thus increasing the temperature of the hair follicle. When the temperature rises to a certain level, the hair follicle will be irreversibly damaged, and the hair will lose its original environment and be completely removed.

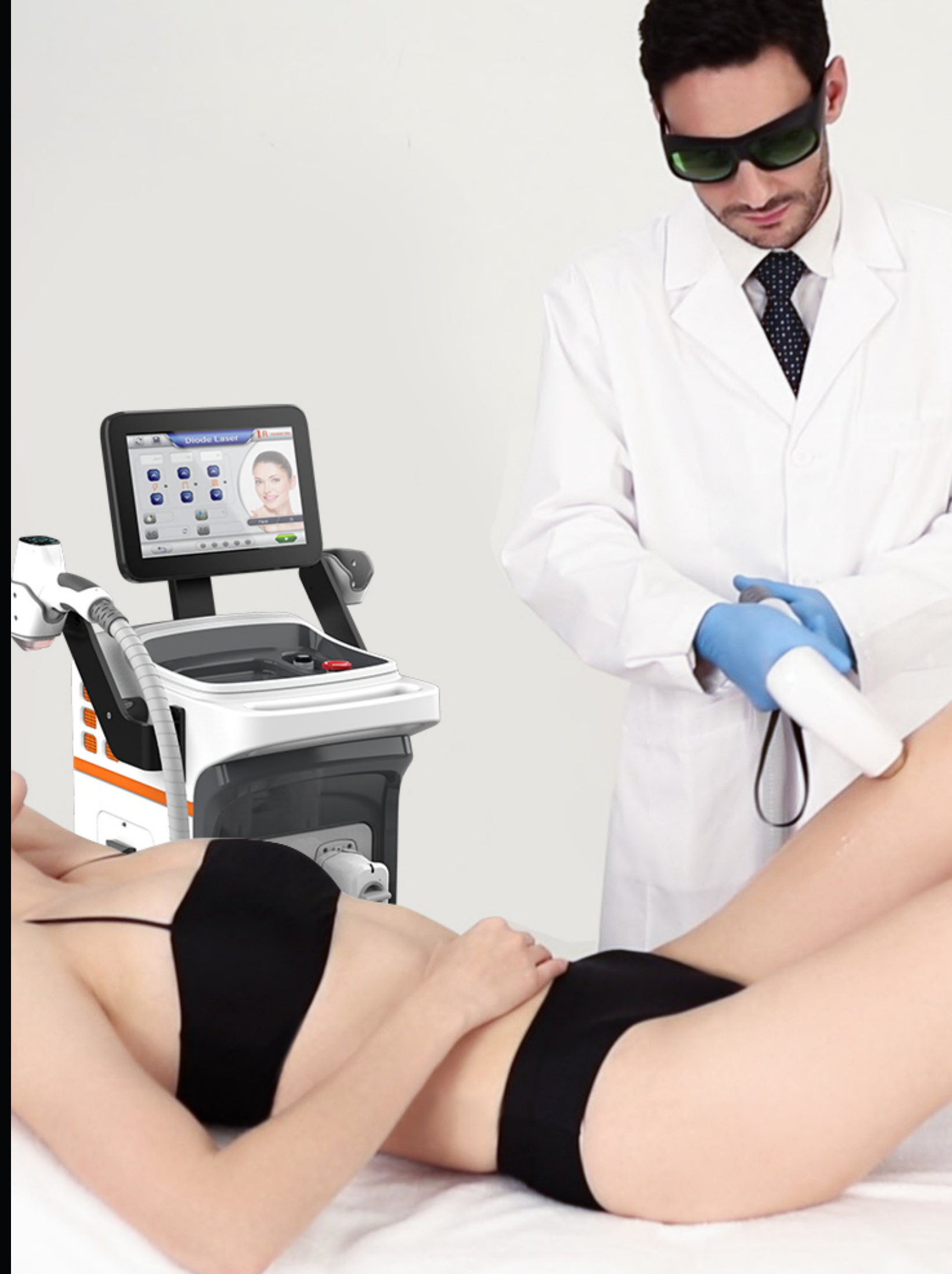

FOUR WAVELENGTHS, BEST RESULT FOR ALL SKIN TYPES AND ALL HAIR TYPES

755nm

755nm offers more powerful energy absorption by the melanin chromophore, making it ideal for light colored and thin hair. It targets the Bulge of the hair follicle and is especially effective for superficially embedded.

808nm

808nm is classic wavelength in laser hair removal, the 808nm wavelength, offers deep penetration of the hair follicle with high average power, suitable for most skin types and hair types.

940nm


This wavelength provides better absorption by oxyhaemoglobin (coagulation of microvessels that feed the hair root) and a moderate penetration depth that makes it suitable for treating all hair types. Melanin absorption is low, making it very safe for working with dark phototypes.

1064nm

1064nm wavelength is characterized by lower melanin absorption, making it a focused solution for darker skin types. At the same time, the 1064nm offers the deepest penetration of the hair follicle.

755nm
808nm
940nm
1064nm



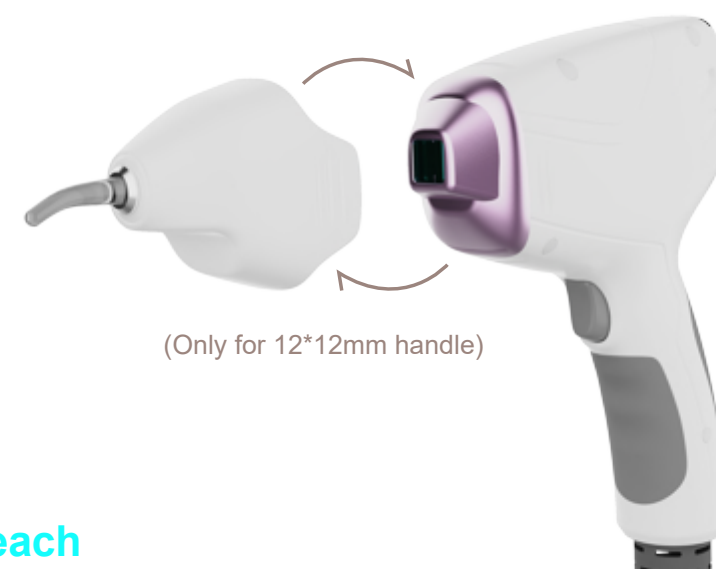

FACIAL TIP

Reaching the hard to reach

The special Facial TIP allows treatment of hard to reach areas, including the ears, nostrils and glabella.

AMERICAN COHERENT LASER BAR



100 million shots warranty

Crystal and laser are packaged in one, with low energy loss and long life

TOUCH SCREEN HANDLE

More than just display

1.54-inch high-definition OLED touch screen, treatment parameters can be adjusted directly on the handle screen.



Energy adjustment




Pulse width adjustment




Frequency adjustment



Real-time temperature monitoring

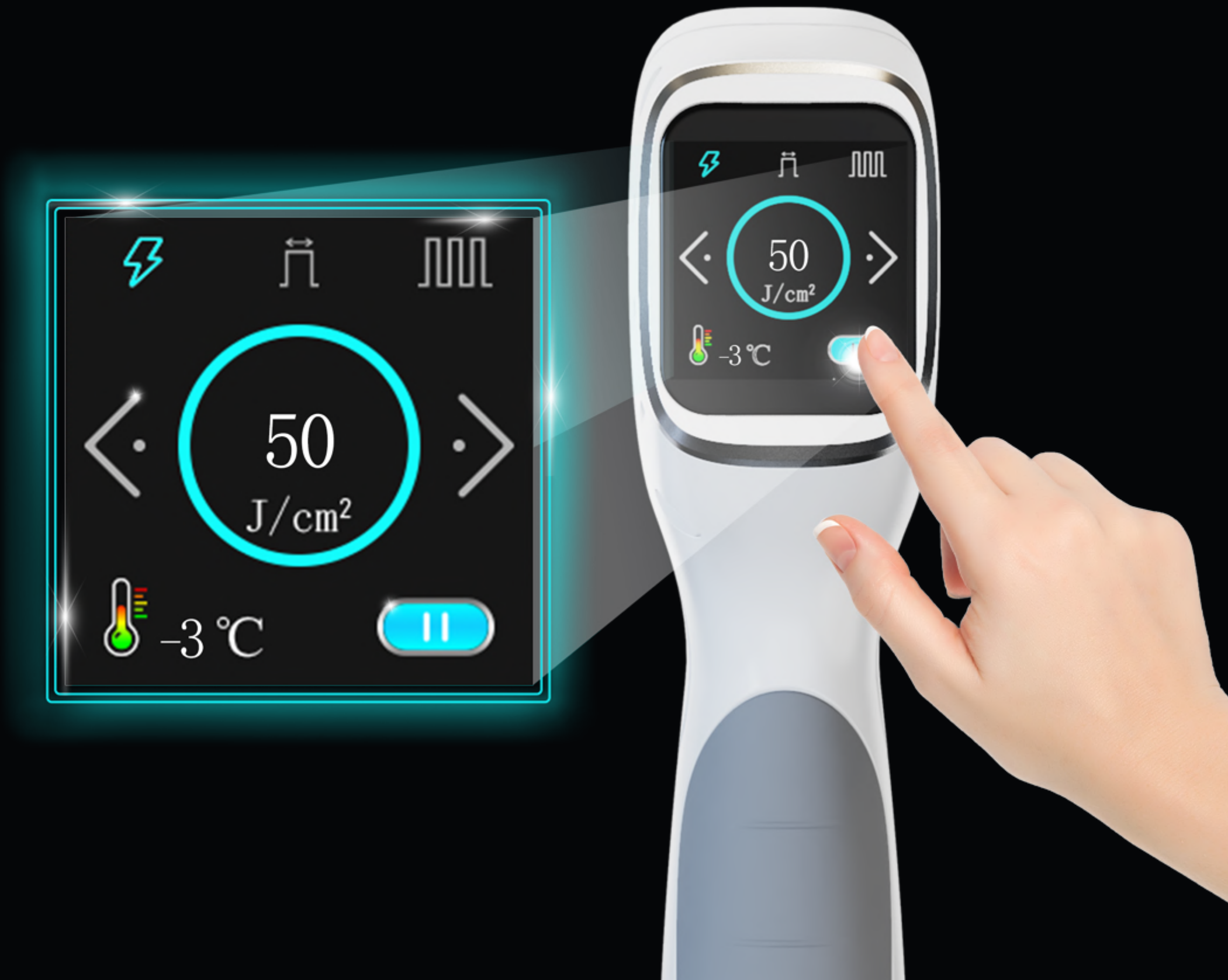

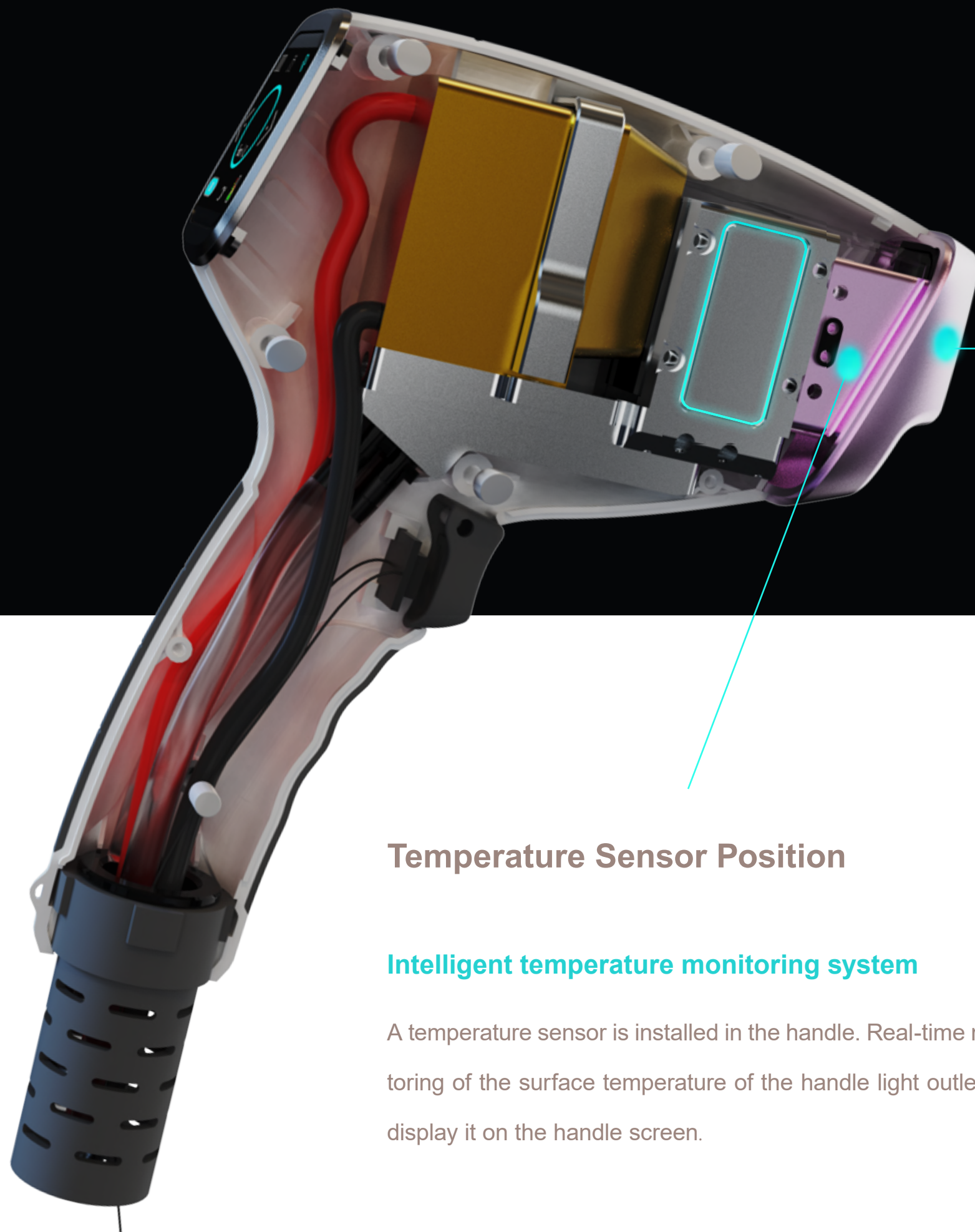

Sapphire crystal

Best light and cold conductivity, ensuring the effectiveness and comfort of the treatment



Super refrigeration system

Two semiconductor cooling plates, 108w super cooling power, the surface temperature of the spot can reach -30 °C

Temperature Sensor Position

Intelligent temperature monitoring system

A temperature sensor is installed in the handle. Real-time monitoring of the surface temperature of the handle light outlet and display it on the handle screen.



INTELLIGENT OPERATING SYSTEM

The operator only needs to select the part and skin color to be treated, and the system will automatically generate the recommended treatment parameters. Even an inexperienced operator can easily grasp the operation of the instrument.

7 Different Languages Interface Available




20 REAL-TIME INTELLIGENT MONITORING

The whole machine is equipped with 20 intelligent detection functions, including water flow detection, water temperature detection, handle temperature detection, handle connection status detection, water level detection, cooling sheet work detection, etc.

Handle connection status detection

Cooling sheet work detection, etc

Water level detection

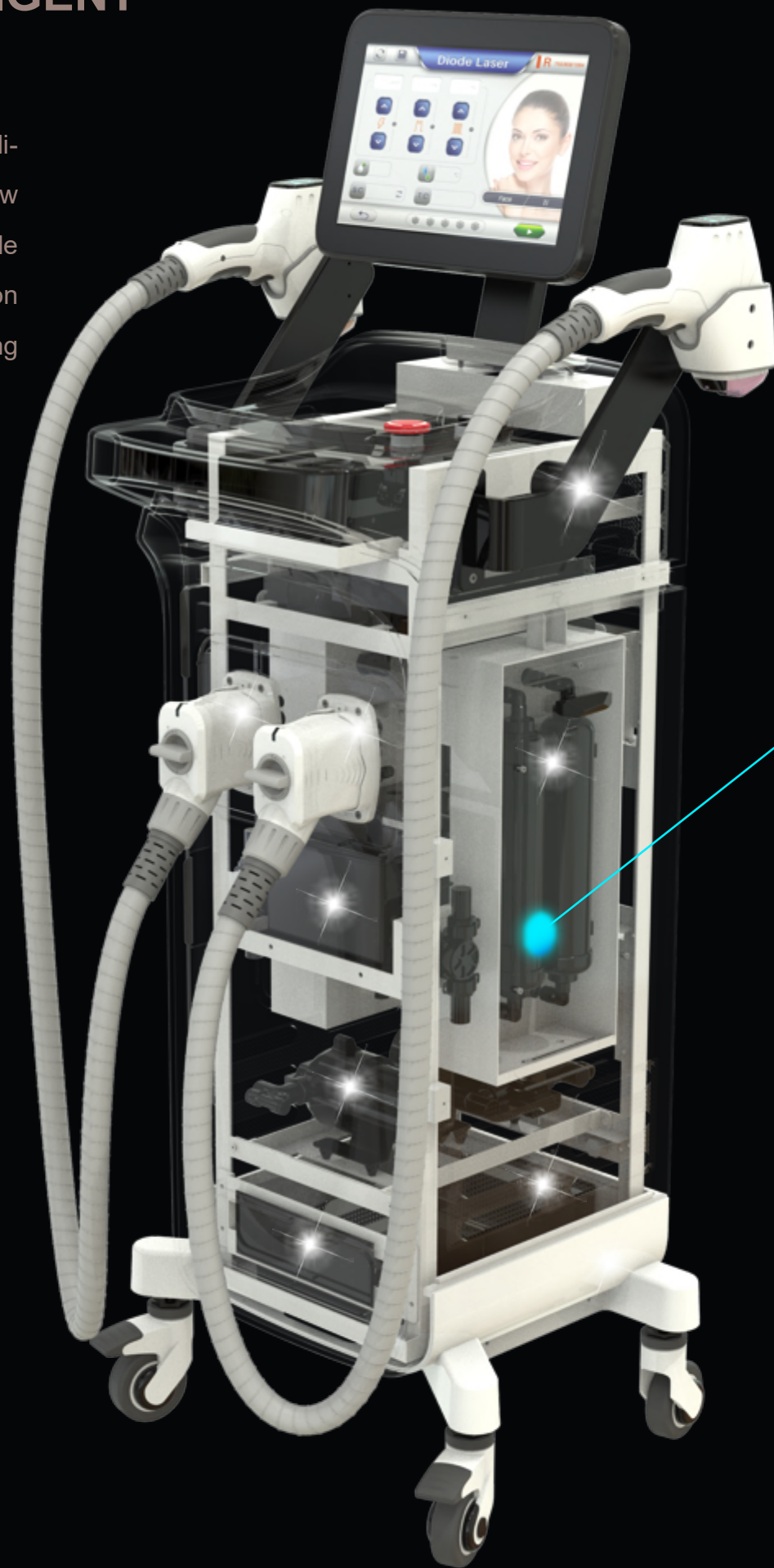

Water flow detection

Handle temperature detection

Water temperature detection

Water level detection

Particulates
0.01-0.1µm

Metal Ion

PP filter

Ion filter

Particulates
0.01-0.1µm

Metal Ion

PP filter

Ion filter

4 FILTERS, QUADRUPLE PROTECTION

Each handle is fitted with an advanced two-stage water filtration system. The first stage adopts PP cotton to filter out impurities and prevent laser blockages, the second stage uses a special Ion filter to filter out metal ions, avoiding inner laser corrosion and prolonging system life.

DIOLASHeer*ice*



SPECIFICATIONS

Power	3000W
Power of handle	1200W+600W
wavelength	755nm+808nm+940nm+1064nm
Machine screen	12.1 inch
Handle screen	1.54 inch
Energy density	1-120J/cm2 （Deviation $\leq\pm 2\%$ ）
Pulse width range	1-200ms
Spot size	12*12mm ; 12*20mm;
Frequency	1-10HZ
Cooling System	TEC cooling system
Net Weight	57kg
Dimension	470*500* 1330mm
Package Size	530*492*1120mm
Fuse specification	Ø5×25 10A
Voltage	AC220V $\pm 10\%$ 10A 50HZ , 110v $\pm 10\%$ 10A 60HZ

