

\$20,500
Includes Free Shipping

NEW GENERATION
SUPER DIODE LASER
HAIR REMOVAL MACHINE

Multi

Multi Spot Sizes

2000w

2000w High Power

Short


Super Short Pulse Duration

DIOLASHeer *ice* Multi



INTELLIGENT TEMPERATURE MONITORING SYSTEM

Real-time monitoring of the surface temperature of the handle.

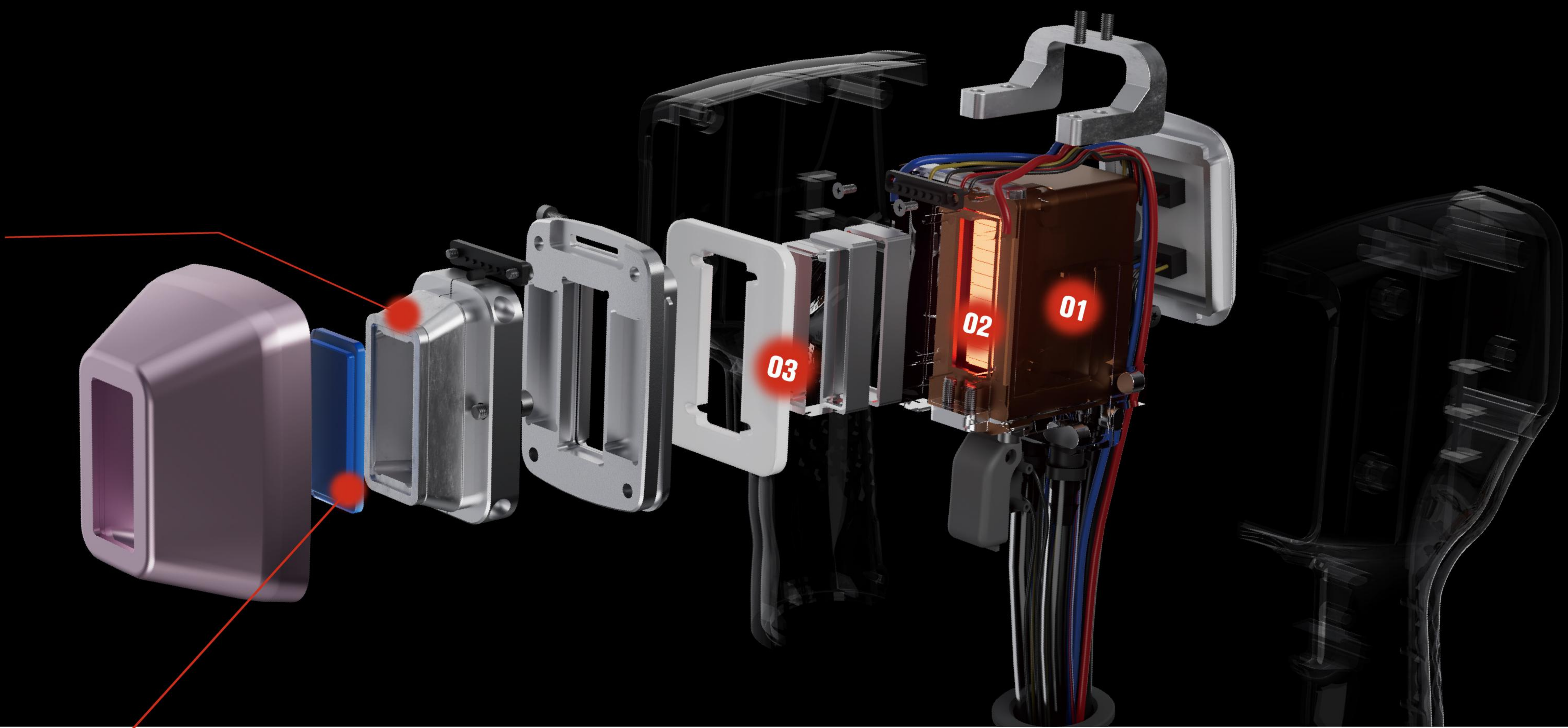

SAPPHIRE CRYSTAL

Best light and cold conductivity, ensuring the effectiveness and comfort of the treatment.

01



AMERICAN COHERENT LASER BAR

100 million shots warranty.

03



360° COOLING PELTIER

108w super cooling power, the surface temperature of the spot can reach -30°C .

02



COLLIMATORLENS FOCUS BEAM

Collimate the light path, reduce the divergence angle of the beam, and turn the divergent light path into a parallel light path.

Four Wavelengths Optional

FOR ALL HAIR TYPES AND SKIN TYPES

755nm

For the fine and residual hair.



808nm

For general laser hair removal at a super speed.



940nm


Effectively penetrates deep and absorbed by the target chromophore.



1064nm

Deeper follicle penetration and darker skin focused.



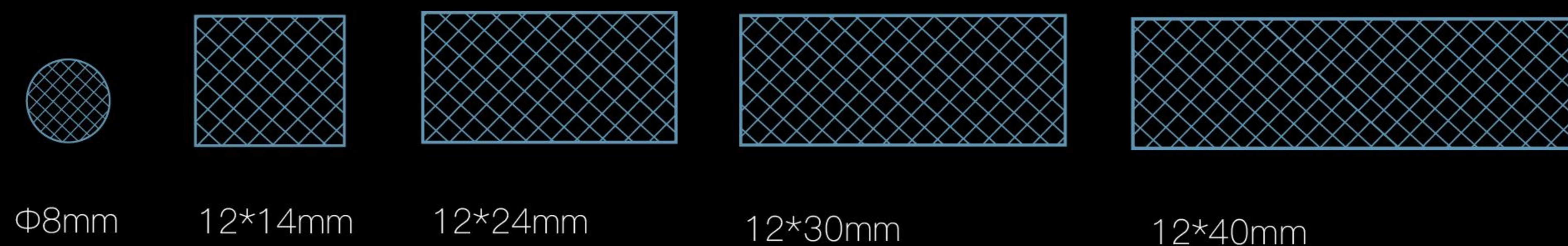

2000W HIGH POWER



MULTIPLE SPECIFICATIONS TREATMENT HEAD

A variety of different powers and different spot sizes are available to meet the treatment needs of different parts of the customer.



TOUCH SCREEN HANDLE


1.54-inch high-definition OLED touch screen, treatment parameters can be adjusted directly on the handle screen.



Frequency Adjustment




Energy Adjustment

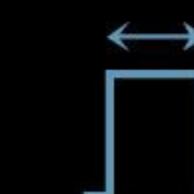

Pulse Width Adjustment

Advantage

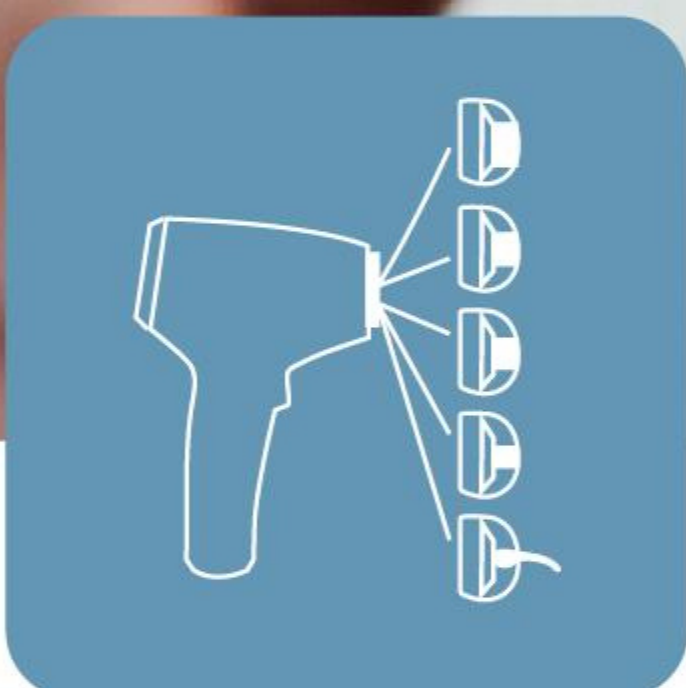

IDENTIFIABILITY

The system can recognize different treatment heads automatically.

750G


LIGHTWEIGHT

The handle weighs only 750g, with humanized design, light and labor-saving





360° COOLING PELTIER

Cooling reaches -27.0°C in 2 minutes.

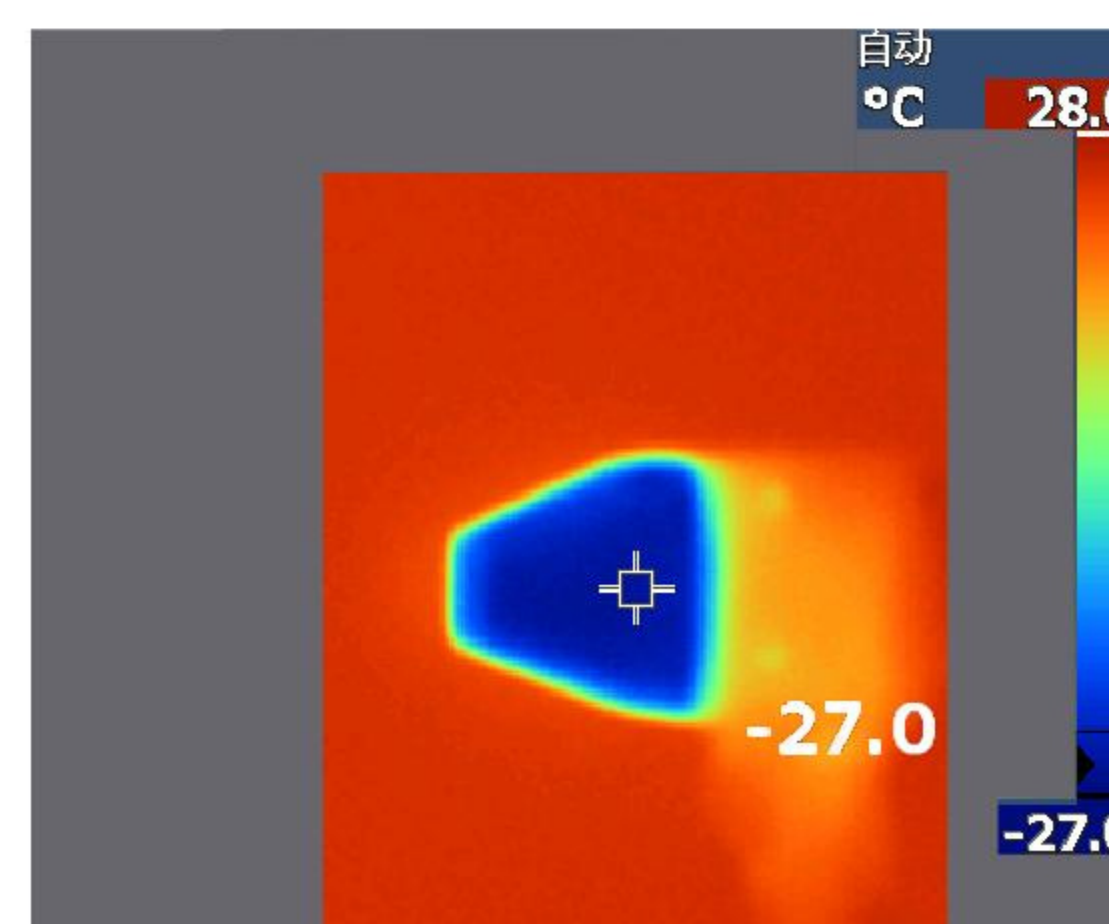

UNIFORM LIGHT SPOT

Uniform light spot, more safer and effective.

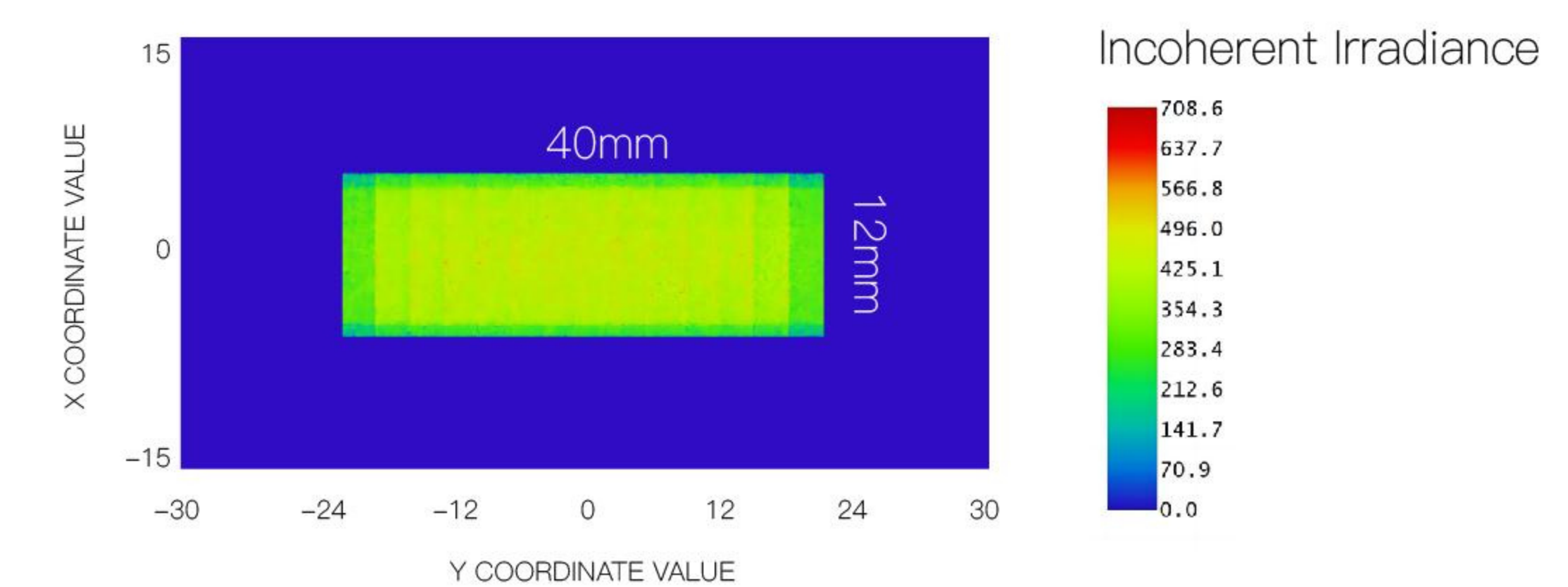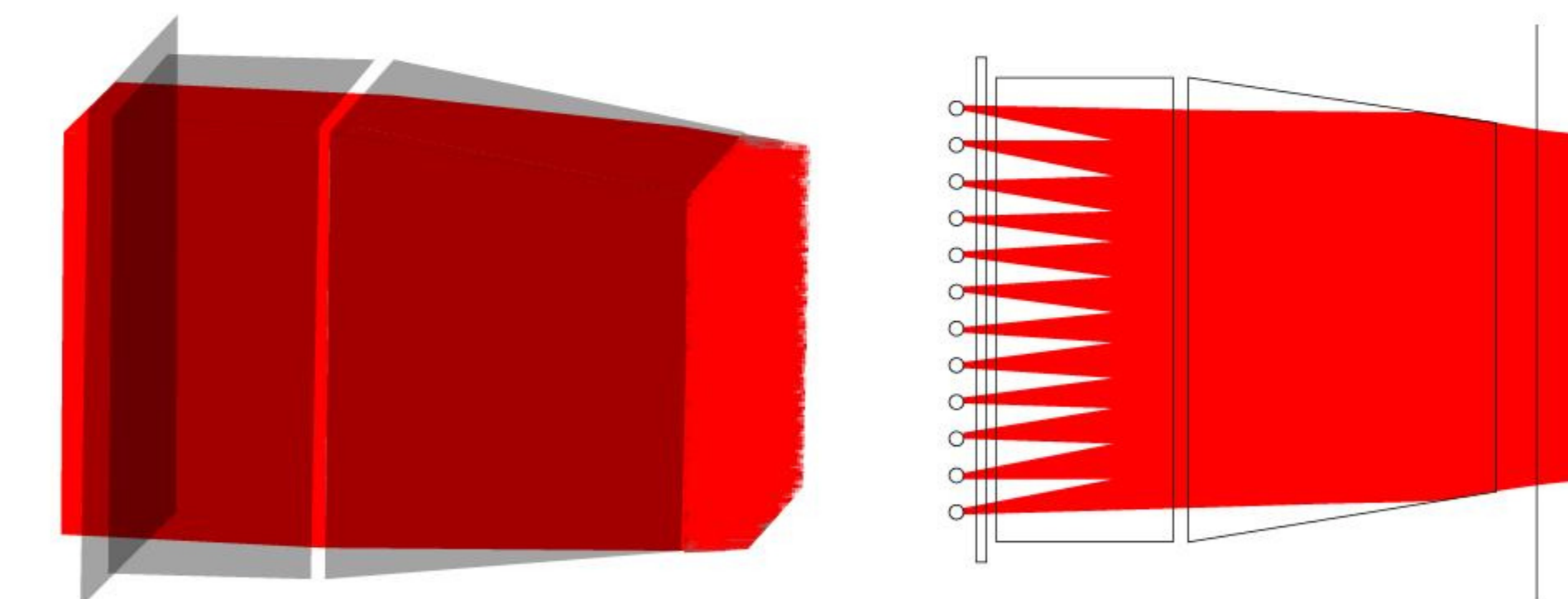

HAIR DETECTION TO INTELLIGENTLY



Unique al diode laser technology, equipped with hair detection the results of each treatment are visible to be seen.



Before



After one treatment



After two treatments



More Effective



Safer



More Professional

ACTUAL EFFECT OF DEPILATION

LEGS



Before




After

ARMPIT



Before

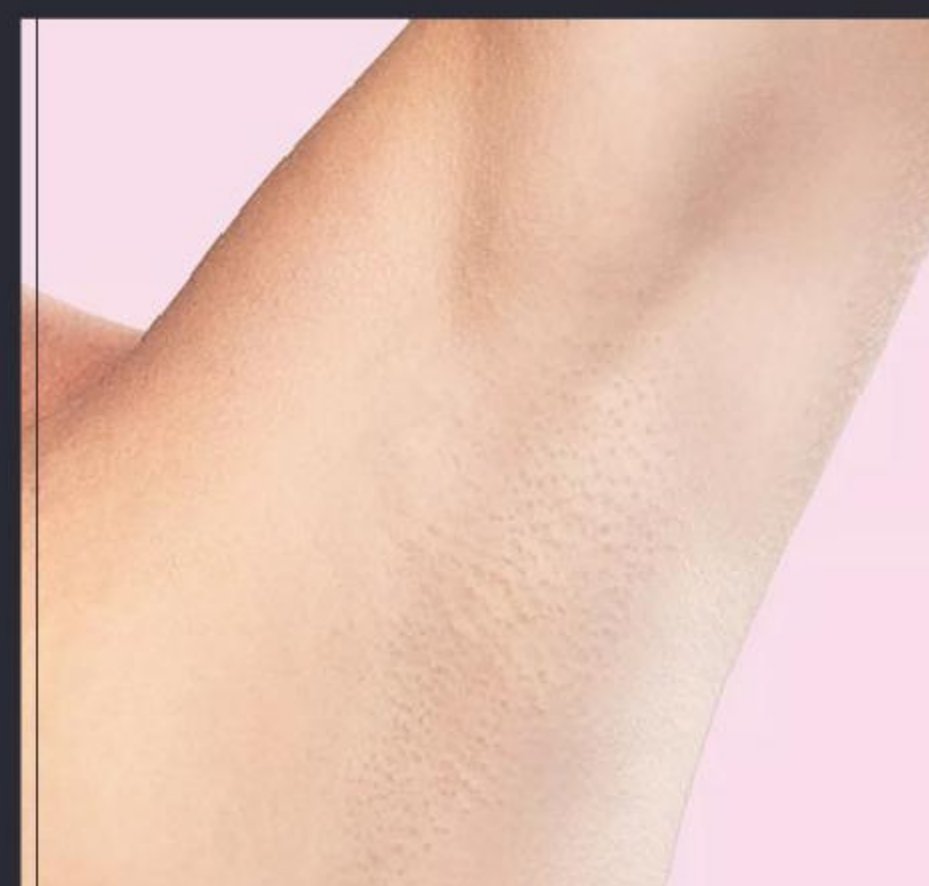

After

BACK



Before




After

BIKINI




Before

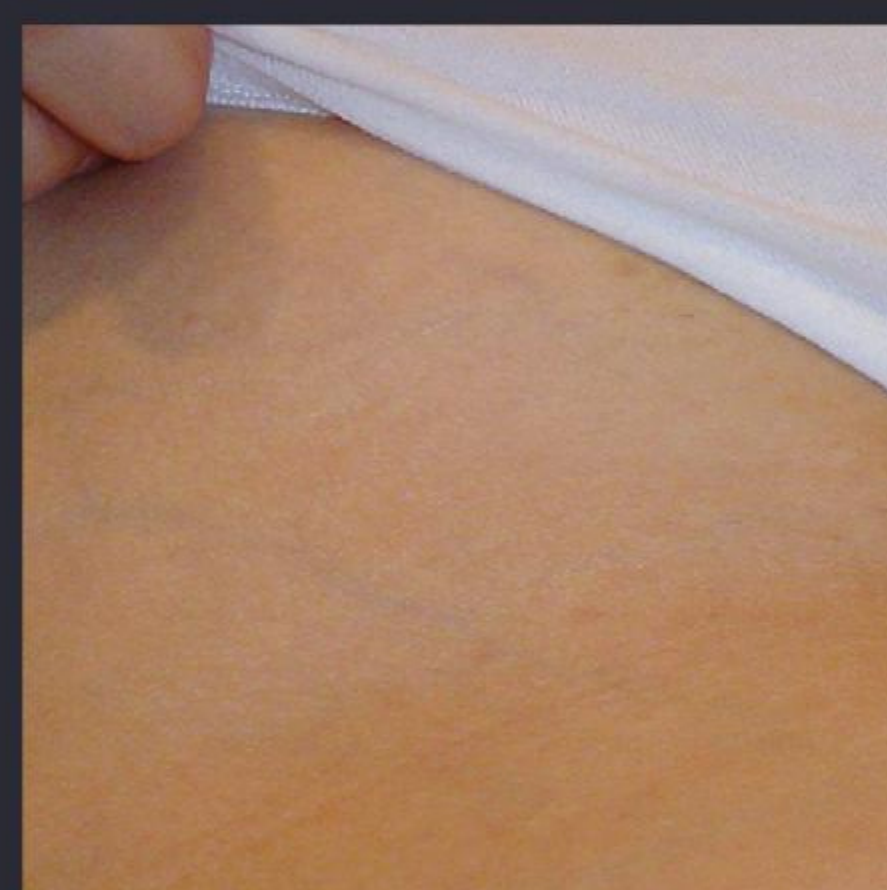


After

CHEEK

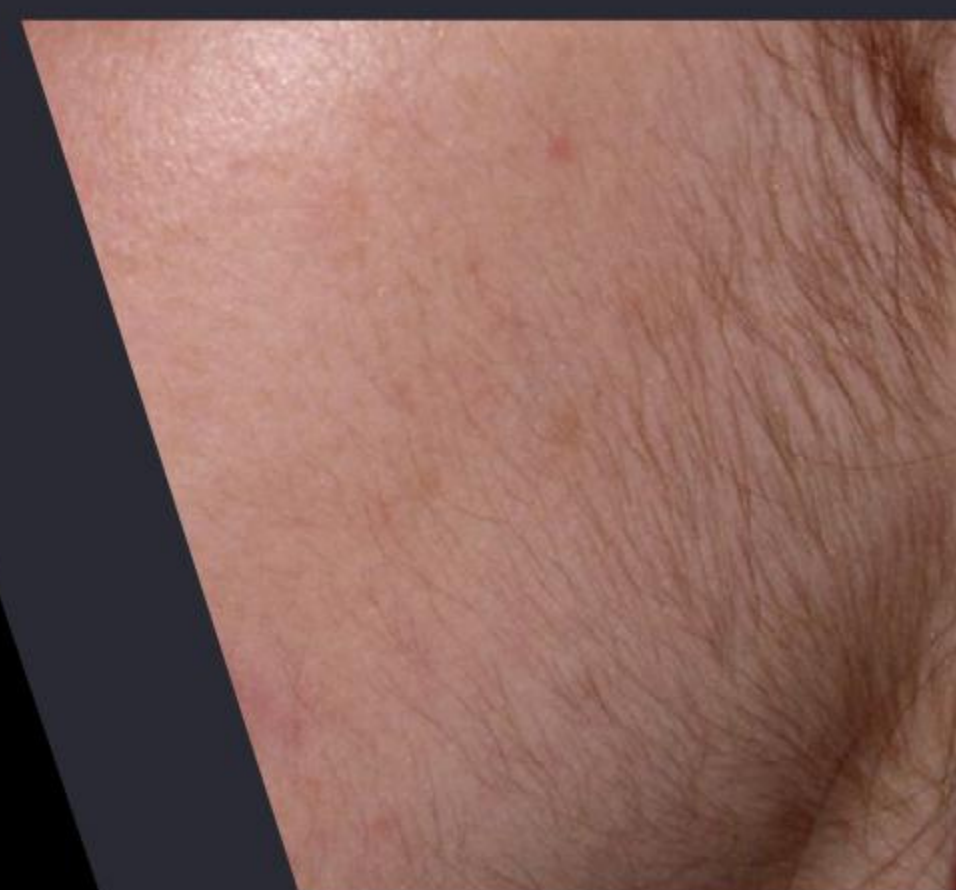


Before

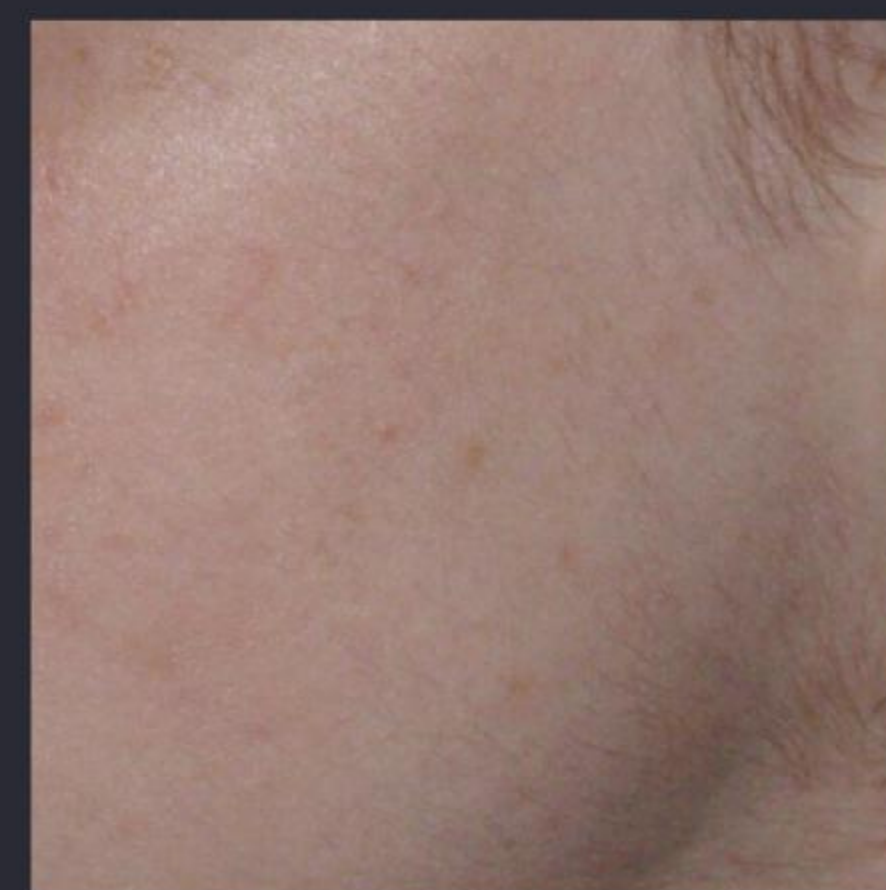


After

LIP

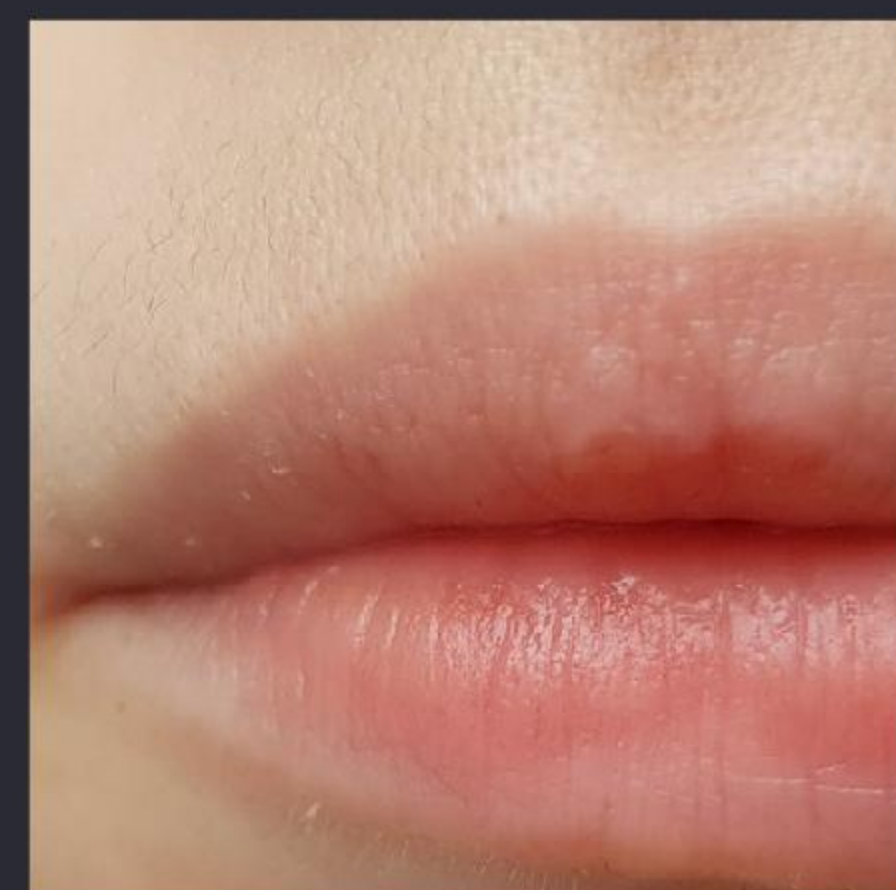

Before

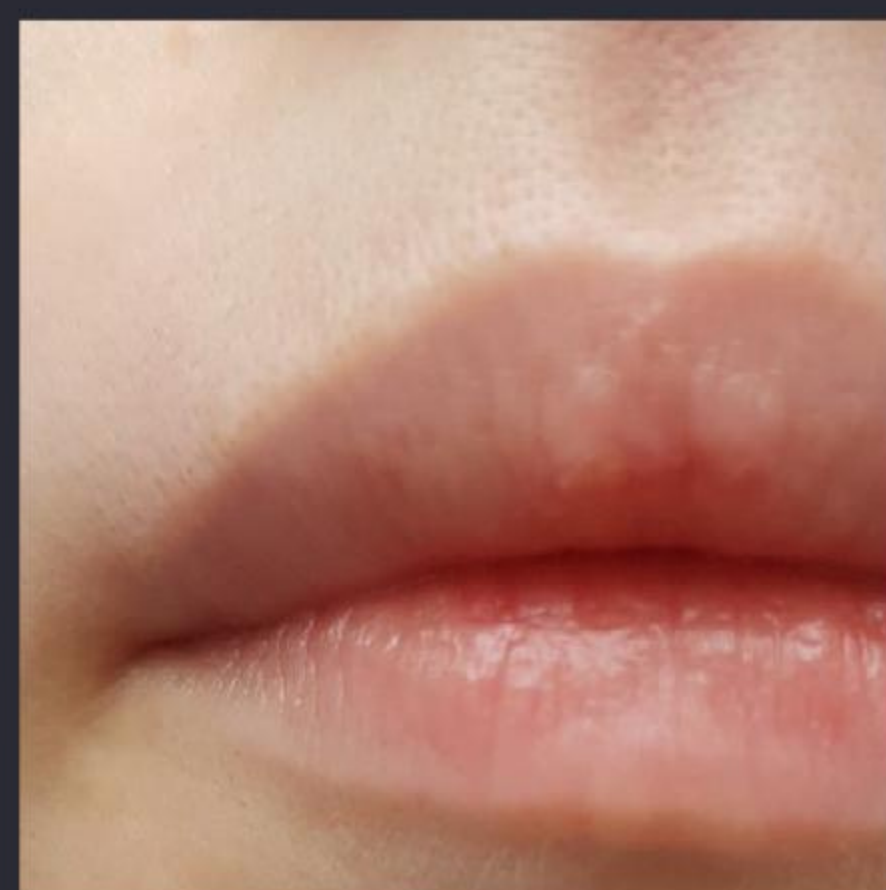

After



ANDROID SYSTEM

WIFI AND BLUETOOTH CONNECTION

CLOUD DATA TRANSFER

AUTOMATICALLY GENERATE TREATMENT PARAMETERS

After selecting gender, treatment area, skin color, hair color, and hair density, the system will automatically generate treatment recommendation parameters.






REMOTE CONTROL SYSTEM

Remotely Control Machine On Mobile

Remotely Setting Parameters On Mobile

Remotely Create And View User Information On Mobile

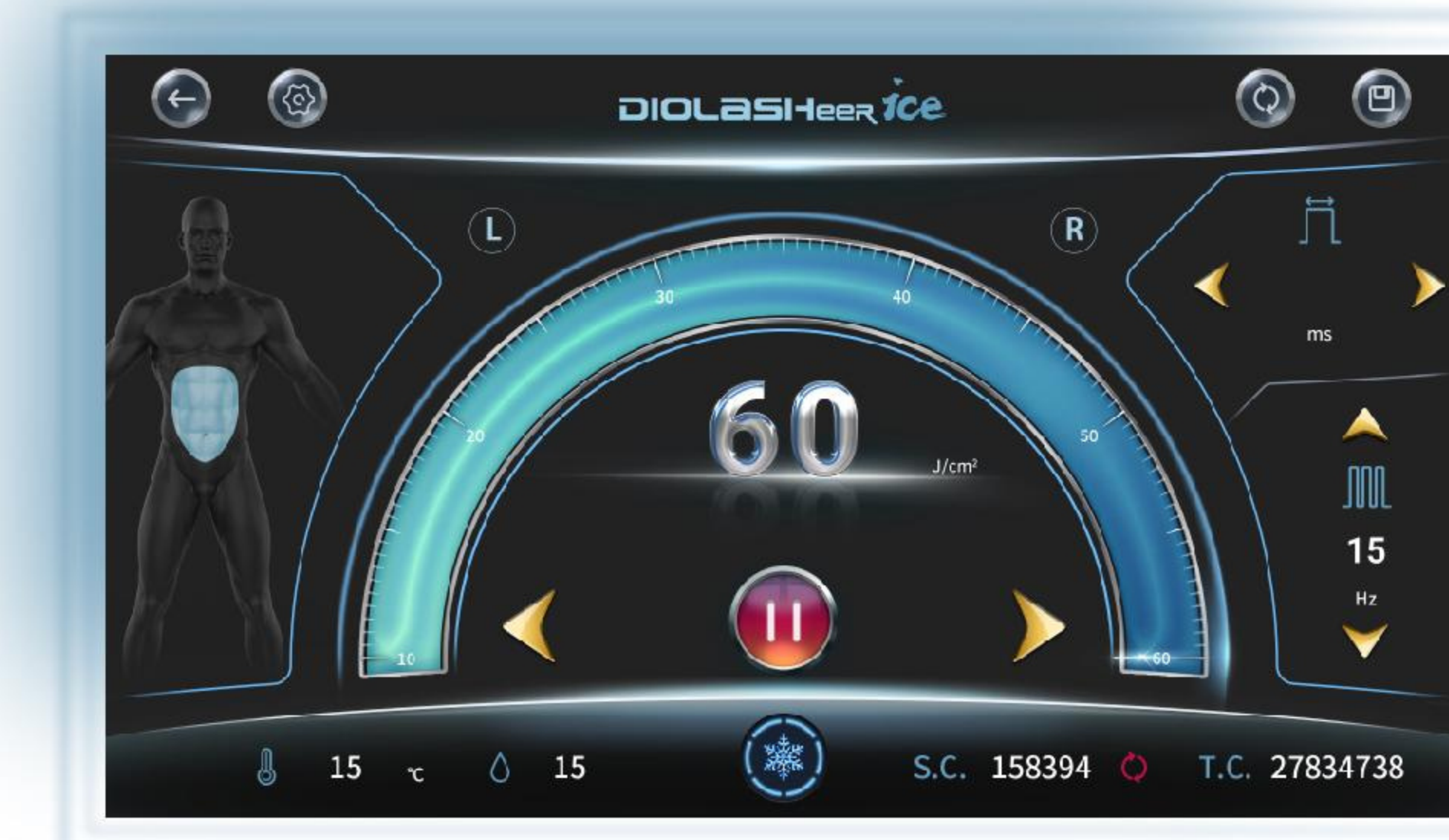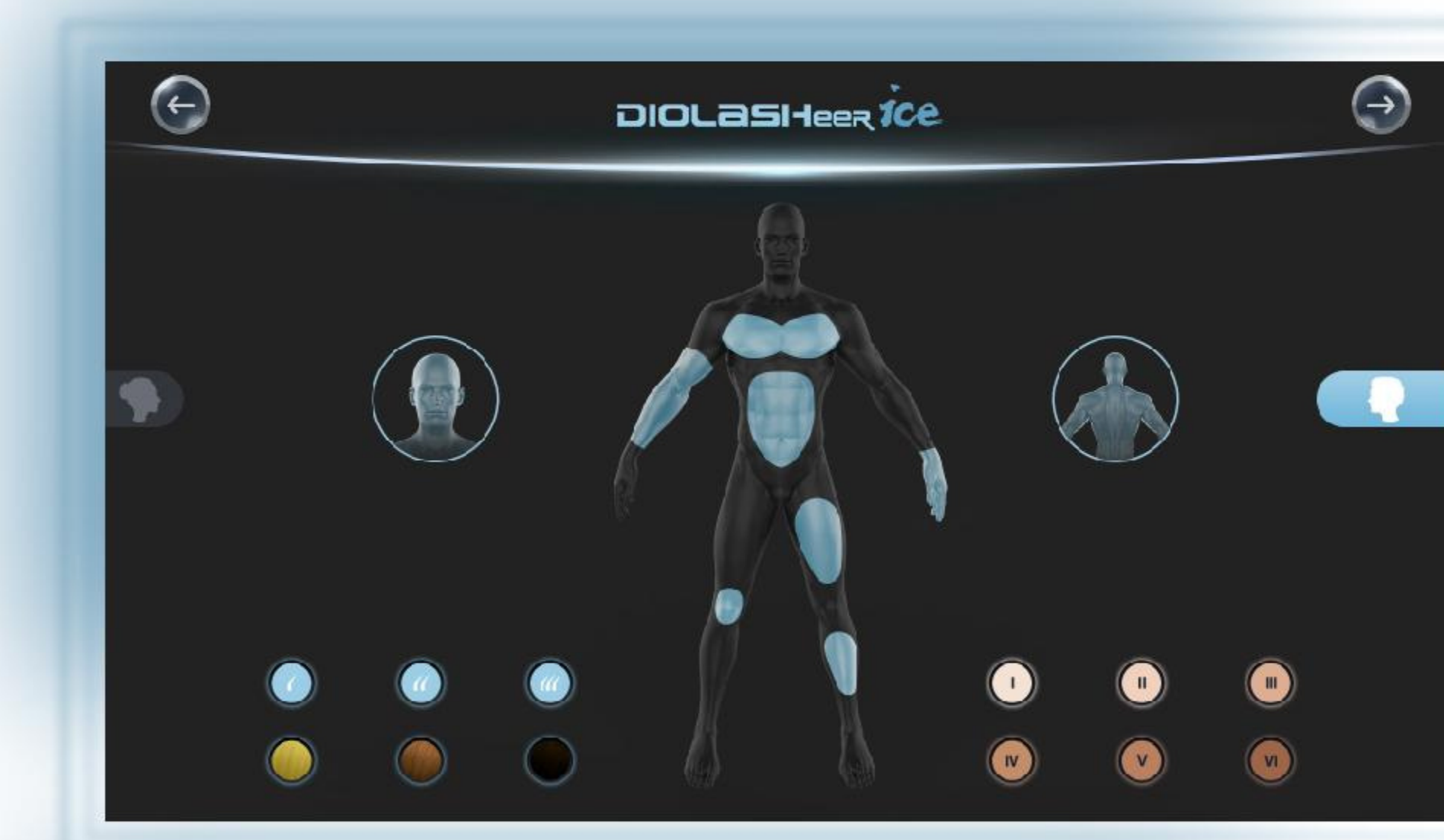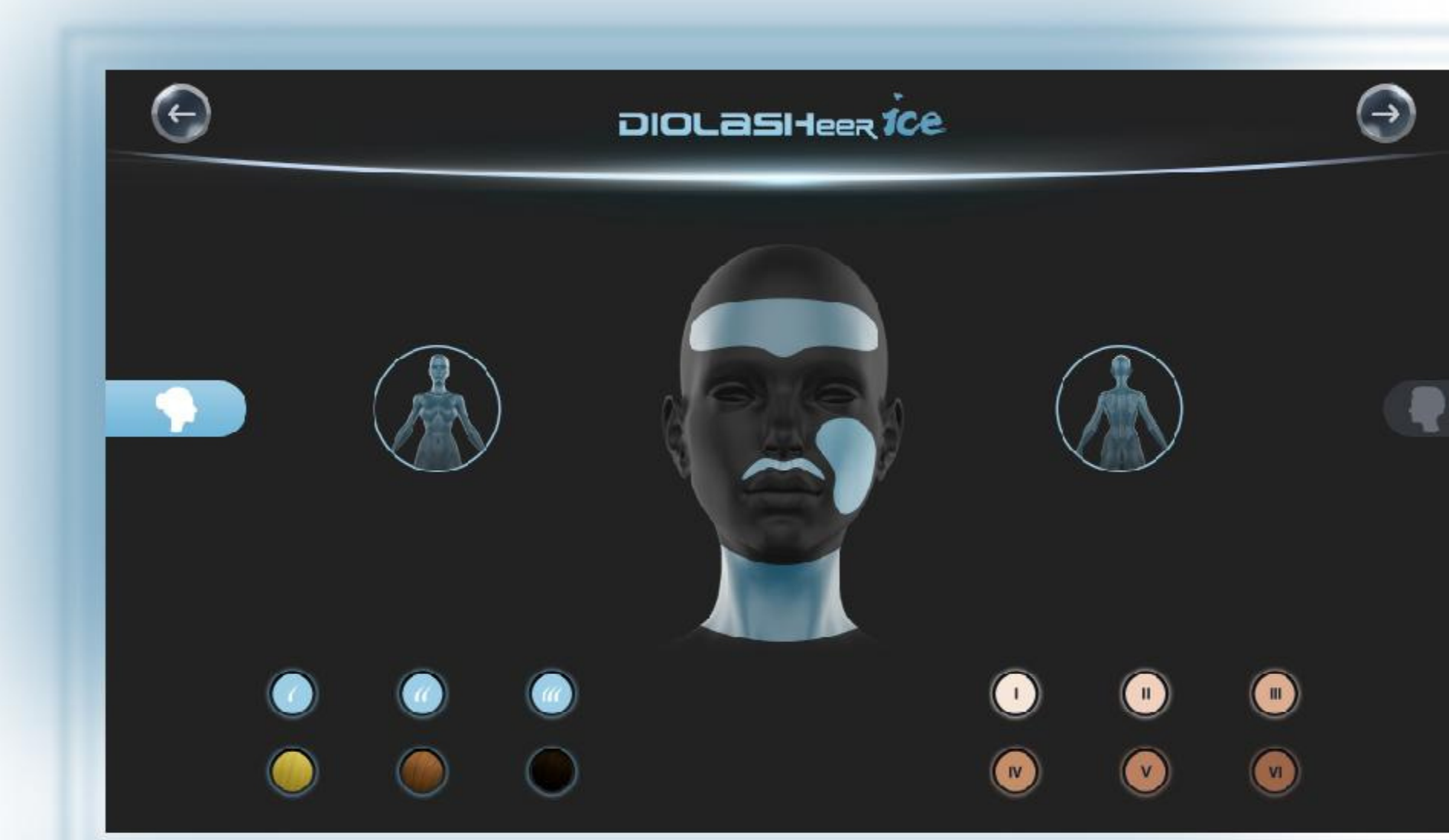

SPECIFICATIONS

Power	3000w
Output Power	2000w
Wavelength	808nm 755nm 940nm 1064nm Optional
Energy Density	10–60 J/cm2 (Deviation≤±2%)
Pulse Width Range	1–200ms
Spot Size	Φ8mm 12*14mm 12*40mm Standard; 12*24mm 12*30mm Optional
Frequency	1–20hz
Cooling System	Compressor Cooling System
Voltage	AC220V±10% 10A 50Hz , 110V±10% 10A 60Hz